R R Hyperbaric Oxygen Chamber



Model:ST801

Chamber material

Abrasion and pressure resistance **Good resilience**

Double thickened window

Anti-rotation base

SN/PD sticker

Fixing the chamber to prevent rotatation

MACY-PAN® 801

hyperbaric O₂ chamber

Transparent window

Easy to see wide field of view



Emergency pressure relief valve

Oxygen inlet

Spare valve

Air inlet

Internal pressure gauge

Real-time detection of chamber pressure

Manual pressure reducing gauge

Operated both internally and externally

External pressure gauge

Double-length internal and external zipper



Automatic constant pressure valve

Product parameters



Chamber size	80*225cm(31*88inch)		
Set pressure	1.3ATA/1.4ATA/1.5ATA		
weight of chamber	about 13KG		
hamber material	TPU		



Oxygen concentrator (Supply Oxygen)



Air Cooler

(Cooling and dehumidification)





Hyperbaric Oxygen Chamber

Air compressor



The product offers sufficient power, ensuring efficient operation. It incorporates a built-in filter cotton that effectively purifies the compressed air, removing impurities and contaminants. This helps to maintain the quality of the air output, improving the performance and longevity of the product. Furthermore, the product is equipped with a silencer that effectively reduces noise levels. These specialized features enhance the overall functionality and user experience of the product.

Air conditioner ((optional))

Air conditioning can maintain a suitable temperature inside the hyperbaric oxygen chamber, providing a comfortable environment for people.



Air cooler



Low-cost cooling and humidity reduction equipment for a more comfortable working environment.


Imported high performance molecular sieve effectively filters 99% of bacteria, providing you with ultra-pure forest-like oxygen.



Oxygen Concentrator



All-in one machine (optional)

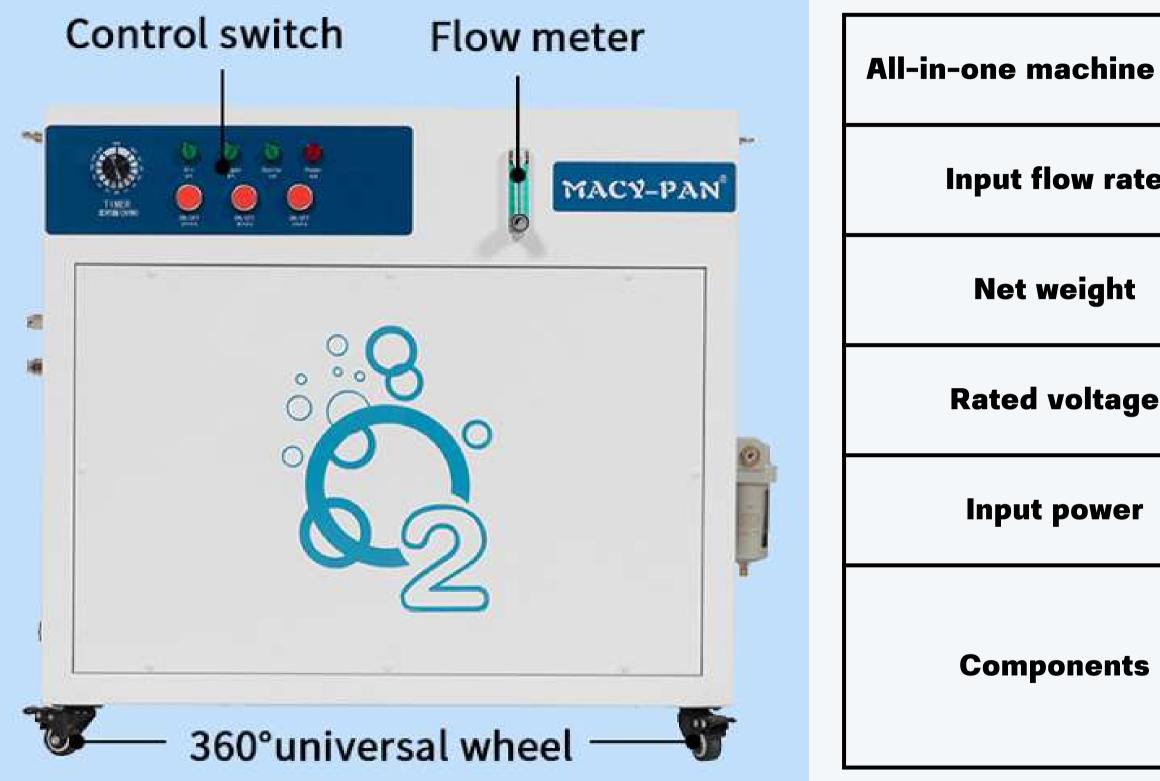

e size	76*42*72cm	
e	70L/min	
	about 80KG	
Ð	220V/110V	
	1280W	
	Air compressor, Oxygen concentrator, Air cooler	
		Y

Technical core





Manual pressure reducing valve

The manual pressure reducing valve can be adjusted both inside and outside.

Automatic exhaust Emergency pressure relief valve

In an emergency, you can press th emergency pressure relief valv quickly decompress.





he		Automatic	exhaust	when	the		
ve	to	pressure	reaches	the	set		
		pressure to keep the chamber					
		pressure constant.					




TPU chamber material

Nano composite fiber. Maintain very good stability in long time or high intensity use environment.



Dual pressure monitoring system

Internal and external pressure gauges monitor real-time pressure and can also be observed by the intelligent control display panel.

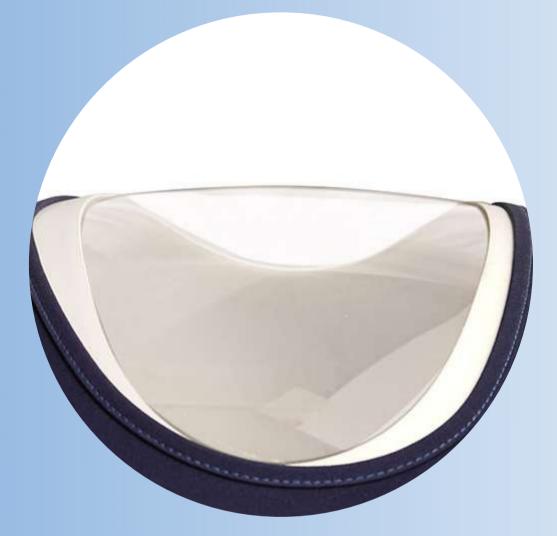

Transparent window

Can see the situation of the people inside the chamber. Full light to eliminate tension.

Multiple oxygen intake system



1. Through oxygen headset (about 93% oxygen)





Through oxygen 2. mask (about 93% oxygen)

4. Breathe oxygen directly from the chamber without using any instruments (oxygen concentration in the chamber is about 30% to 40%)



Through Nasal oxygen 3. tube(about 93% oxygen)

Packing and shipping



Our customers







Customer feedback



27分钟前

展开

(7)

Ŝ.

Dear Sherry,

r

Yesterday, I received the hyperbaric chamber! Everything went very well with the shipping and customs. When the packages arrived, I was amazed at how well and carefully everything was packed! I give you a 5 star rating (the maximum) for the shipping and packaging! When I opened the boxes, I was so happy to discover the excellent guality of your products!!!! I checked everything! The materials that you use are of very good quality. You are really professionals!!!! Please forward this email to your CEO. My sincere congratulations to such a company! Five stars for the product quality!!!! The maximum! But the most important is the magnificent customer support that you provide! I am telling you the truth. 1 billion stars for such a service!!!! Congratulations on such outstanding customer service. Because of all these, please be sure, that I will recommend you to all my friends!!!

Thank you so much for showing me through videos and PDFs how to install the chamber and how to operate it. I am really, really happy with my purchase!

Thank you also for your kindness in promising to assist me further in the future.

All my best wishes from Sibiu, and I wait for all of you to come to visit Romania!!!

> の回复 ⇔转发

MACY-PAN®

ST801 Lying Type Portable Hyperbaric Oxygen Chamber



client feedback

America customer World University Skiing Champion Natalie Good She uses chamber for exercise recovery.

Sarah Milne ireat to have Natalie Good in the hyperbanic chamber . Nats is the fire and only New Zealander to have won a gold medal at the World Universit Games, when she won the Slope Style at Trentino in Italy. Nats is just back from her uber successful trip to compete in Euro







MACY-PAN®

ST801 Lying Type Portable Hyperbaric Oxygen Chamber



New Zealand customer famous ultramarathon expert Lisa She uses chamber for her mum's recovery after Aneurysm surgery. Mum recover very fast and now can walk by herself.







五0 円0

Certificat

證明書

0 O

Certific

